


PICC-Nate[®]

Peripherally Inserted Central Catheter
For Neonatal Patients

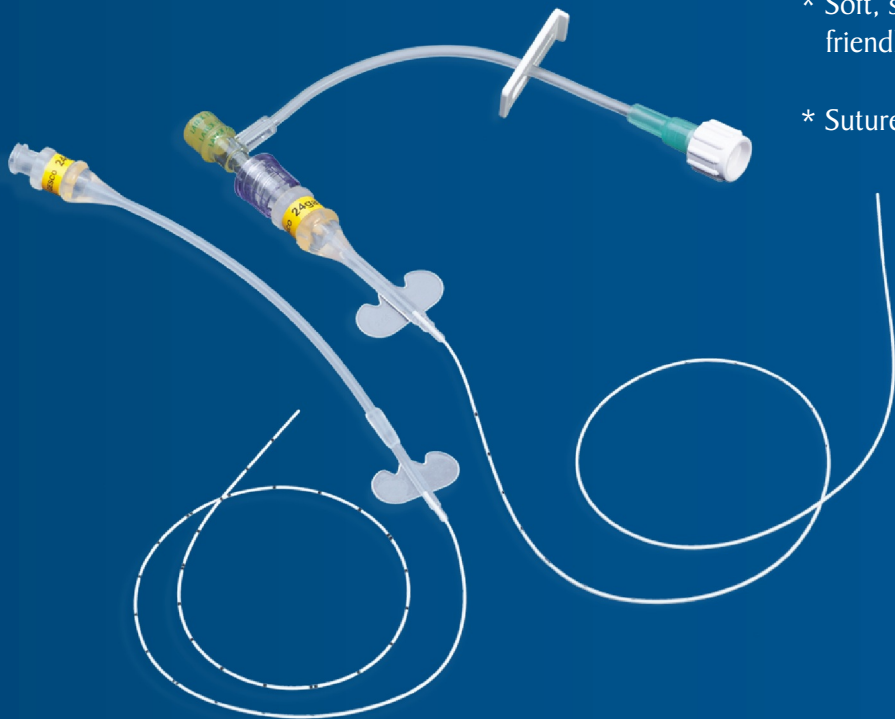

- * Slender, low-profile hub reduces patient discomfort.
- * Soft, strong silicone shaft material provides long-term friendly care.
- * Suture wing is profiled for excellent hub securement.
- * Catheter flow rates provide excellent infusion therapy.
- * Accurate catheter positioning accomplished with depth marking at every centimeter.
- * Reduced catheter shaft diameters minimize puncture wound.
- * Barium sulfate-loaded catheter body provides radiopacity for fluoroscopic visualization.

Specifically designed to minimize trauma to the neonatal patient.



- * Double-beveled needle tips allow precise vessel entry.
- * Funneled entry guides PICC introduction.

Utah Medical Products, Inc.
provides a broad range of critical
care devices including the Gesco[™]
family of neonatal products.



All Gesco products are Latex-free.

Product Description

PICC-Nate catheters and introducers are available in a variety of kit configurations, in "Basic" or "Complete" pre-assembled procedure trays.

Basic Catheter Kit Contents: Catheter and Introducer only.

Complete Catheter Tray Contents:

| | |
|--|--|
| 1 Single Lumen PICC Catheter | 1 each: Scissors |
| 1 Peel-Away Introducer | 1 each: Syringe with 18 Ga Needle |
| 2 each: Tape measure (1/8" x 15") | 1 pkg: Alcohol Swabsticks (3 per pack) |
| 1 pkg: Povidone Iodine Swabsticks (3 per pack) | 1 each: Extension Set with Injection Port (excluding PC-2LPS, PC-3LPS) |
| 8 each: 2"x2" Gauze Pads | 3 each: Tape Strips |
| 8 each: 4"x4" Gauze Pads | 1 each: Instructions For Use |
| 1 Fenestrated Drape with Tape | 1 each: Poly-Lined Drapes |
| 1 each: Tourniquet | 1 each: Micro Forceps (Half Curved Fine Point) |
| 2 each: Point Lok Devices | |

Ordering Information and Specifications:

Complete Catheter Trays

| Catalog Number | Box Quantity | INTRODUCER | | | Fr. Size | Gauge | CATHETER | | |
|----------------|--------------|------------|--------|-------|----------|-------|-----------|-----------|-------------------|
| | | Peel Away | Length | Gauge | | | Material | Lumen DIA | Gravity Flow Rate |
| PC-2PS | 10 | X | 25mm | 22 | 1.9 | 24 | silicone | .36mm | .53ml/min |
| PC-2LPS | 10 | X | 25mm | 22 | 1.9 | 24 | silicone | .36mm | .53ml/min |
| PC-2PP | 10 | X | 25mm | 22 | 1.9 | 24 | Tecoflex® | .36mm | .53ml/min |
| PC-3PS | 10 | X | 35mm | 19 | 3.0 | 20 | silicone | .51mm | 2.37ml/min |
| PC-3LPS | 10 | X | 35mm | 19 | 3.0 | 20 | silicone | .51mm | 2.37ml/min |
| PC-3PP | 10 | X | 35mm | 19 | 3.0 | 20 | Tecoflex® | .51mm | 2.37ml/min |

Basic Catheter Kits

| Catalog Number | Box Quantity | INTRODUCER | | | Fr. Size | Gauge | CATHETER | | |
|----------------|--------------|------------|--------|-------|----------|-------|-----------|-----------|-------------------|
| | | Peel Away | Length | Gauge | | | Material | Lumen DIA | Gravity Flow Rate |
| PB-2PS | 10 | X | 25mm | 22 | 1.9 | 24 | silicone | .36mm | .53ml/min |
| PB-2LPS | 10 | X | 25mm | 22 | 1.9 | 24 | silicone | .36mm | .53ml/min |
| PB-2PP | 10 | X | 25mm | 22 | 1.9 | 24 | Tecoflex® | .36mm | .53ml/min |
| PB-3PS | 10 | X | 35mm | 19 | 3.0 | 20 | silicone | .51mm | 2.37ml/min |
| PB-3LPS | 10 | X | 35mm | 19 | 3.0 | 20 | silicone | .51mm | 2.37ml/min |
| PB-3PP | 10 | X | 35mm | 19 | 3.0 | 20 | Tecoflex® | .51mm | 2.37ml/min |

Catheter Introducers

| Catalog Number | Box Quantity | Peel Away | Length | Gauge |
|----------------|--------------|-----------|--------|-------|
| PEEL2 | 5 | X | 25mm | 22 |
| PEEL3 | 5 | X | 35mm | 19 |

Catheters

| Catalog Number | Box Quantity | Material | Gauge |
|----------------|--------------|-----------|-------|
| P-2S | 10 | silicone | 24 |
| P-2P | 10 | Tecoflex® | 24 |
| P-2LS | 10 | silicone | 24 |

Catheter Repair Kits

| Catalog Number | Box Quantity | Fr.Size |
|----------------|--------------|---------|
| P-R2 | 10 | 1.9 |
| P-R3 | 10 | 3.0 |

EU REPRESENTATIVE
Utah Medical Products, Ltd.
Athlone Business & Technology Park
Athlone, County Westmeath, Ireland
+353-90-647-3932

UNITED STATES
7043 South 300 West
Midvale, Utah 84047
(801) 566-1200
Toll Free 1-800-533-4984
Fax (801) 566-2062



UTAH MEDICAL
PRODUCTS INC.



www.utahmed.com

www.order.utahmed.com

P/N 58150 REV 101207